



(EN) PRODUCT NAME

HBsAg Hepatitis B surface antigen rapid test kit (serum/plasma/whole blood). A rapid, one-step test for qualitative detection of hepatitis B virus in human serum, plasma and whole blood. For in vitro diagnostic use only.

PURPOSE

The sample mixing colloid-gold monoclonal antibody moves along the membrane to the T-line and forms the T-line when human serum/plasma and whole blood contain HBsAg according to the principle of the sandwich method with double antibody and gold immunochromatographic test, which is a positive result. The unreacted markers move forward continuously to combine with anti-mouse antibodies to form a control line. If the test line does not appear, it is a negative result.

MAIN COMPONENTS

Sample pad, colloidal gold labeled pad, nitrocellulose membrane, absorbent paper and PVC board. Different components of different batches cannot be used at the same time to avoid erroneous results.

STORAGE AND SHELF LIFE

Store as packed in the sealed pouch at 4-30°C, avoid high temperature and direct sunlight in a dry place, valid for 24 months. DO NOT FREEZE! Some protective measures should be taken during hot summer and cold winter to avoid high temperature or freeze-thaw.

SAMPLE REQUIREMENTS

Use fresh whole blood samples. Whole blood samples with anticoagulant can be stored at 2-8°C for 24 hours and should be used immediately without anticoagulant. Frozen whole blood samples should not be used. (If bleeding occurs, it is recommended to use serum samples.)

MATERIALS

1. Test cartridge
2. Lancet
3. Capillary dropper
4. A packet of alcohol swabs
5. Diagnostic solution
6. Leaflet and method of use

Materials needed but not provided :
- Timer






RESEARCH METHOD

Instructions for use should be read in full before testing. Allow the controls

on the test device to equilibrate to room temperature for 30 minutes (20°C-30°C) before testing. Do not open the inner packaging until ready to use, it should be used after one hour if opened.

Tape:

1. Wash your hands with soap and rinse thoroughly with clean warm water.
2. Leave the box at room temperature before opening it.
3. Open the box, remove the test cartridge and place it on a clean and flat surface. Then remove all other accessories from the box and place them in a convenient place near you. Perform the test within one hour, and you will get the best results if the test is performed immediately after opening the foil pouch.
4. Carefully pull and discard the released lancet cap. Use an alcohol swab to clean the fingertip /middle or ring finger/ that will be pricked. Allow to air dry.
5. Press the lancet on the front side where the cap is removed. The nozzle retracts automatically and safely after use. Massage the hand without touching the puncture site, massaging the hand towards the fingertip of the middle or ring finger to be punctured. Holding the hand down, massage the end that is pricked to get a drop of blood. Without squeezing the dropper bulb, place it in contact with the blood. The blood migrates into the capillary dropper through the capillary to the line indicated on the capillary dropper.
6. Using the dropper, vertically transfer two (2) drops (50-60ul) of the whole blood sample onto the sample test portion (S) of the cartridge, avoiding bubbles. Allow the blood to completely disperse in the well, then clean your finger with alcohol and a swab.
7. Unscrew the cap of the reagent bottle and add one (1) drop (40-50ul) of solution to the test area (S) and start the timer.
8. Observe the test results immediately within 15-20 minutes, the result is invalid after 20 minutes.

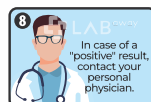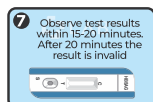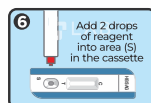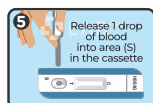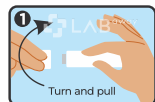

READING THE RESULTS

Positive

Two distinct red lines are displayed. A red line in the control region (C), and another red line in the test region (T). A pink to red line (T), even if very narrow, indicates a positive result.

Negative

No red lines in the control region (C). No lines in the test region (T).

Invalid

No red lines or the control line is not displayed. Insufficient specimen or incorrect procedural techniques are the most likely causes. Review the procedure and repeat the test with a new device. If the problem persists, stop using the kit and contact your local distributor.

LIMITATIONS

The reagent is designed for a qualitative screening test. The concentration of HBsAg cannot be measured by this qualitative test. Reagent results are for clinical reference only, a confirmed diagnosis should only be made by a physician after all clinical and laboratory tests have been reviewed.

PERFORMANCE CHARACTERISTICS

Detection Limit :

The limit should be no more than 2.5ng/mL. Results should be positive when the national HBsAg qualified control samples of adr, adw and ay are detected.

Rejection specificity: Results should be negative when 20 sets of HBsAg negative control samples are detected.

Positive specificity: Results should be negative when 20 sets of HBsAg positive control samples are detected.

The results must be consistent and the degree of staining must be the same when accurate precision control samples from 10 sets of the same batch are detected.

Analytical sensitivity: 1000 mol/L bilirubin, 5.65mmol/L triglycerides, 6.5 g/L hemoglobin has no effect on the result detection. The reagent is not affected by rheumatoid factor, HBsAb, HBeAg, HBcAb. Diagnostic specificity and sensitivity was clinically evaluated comparing the results obtained using the HBsAg hepatitis B surface antigen rapid test kit (whole blood) and EIA HBsAg assays.

The results are as follows:

Method/EIA	Evaluation test		Total
	Positive	Negative	
Positive	377	1	378
Negative	0	622	622
Total	377	623	1000

Relative sensitivity: 100%

Relative specificity: 99.84%

Accuracy: 99.90%

PRECAUTIONS

The result is invalid after 20 minutes. Do not use after expiry date. The strength of the qualification control line does not indicate a problem with the reagent qualification, a test result that is clearly visible demonstrates that the reagent is effective. During testing, dark and light color may appear, even very light mark. If the blood is too thick and the climb rate is less than 4 mm/min, add one drop of the solution to the sample pad to help. Do not use other types of qualified control samples to test the reagent. Components from different batches cannot be interchanged to avoid erroneous results.

MANUFACTURER

Qingdao Hightop Biotech Co., Ltd.
Add: No.369 Hedong Road,
high-tech industrial development zone,
Qingdao, Shandong, 266112, China

REPRESENTATIVE

„Griffon Health Solutions Ltd.
gr. Sofia, ul. "Exarch Joseph" №30
Phone: +359 882 399 001
See instructions for use .
Disposable.

Store between 4-30 °C



www.lab-away.com