

(EN) Intended use

The Lab-Away Urine Alcohol Test is a rapid and highly sensitive method for detecting the presence of alcohol in urine.

This test provides only a preliminary screening. To obtain a confirmed analytical result, a more specific alternative chemical method must be used. Clinical caution and professional judgment should be applied to any test result, especially when preliminary positive results are indicated.

Summary

Two-thirds of adults drink alcohol. The blood alcohol concentration at which a person becomes inadequate varies and depends on the individual. Each individual has specific parameters that influence the level of impairment, such as size, weight, eating habits, and alcohol tolerance. Inappropriate and excessive alcohol consumption can be a contributing factor in many accidents, injuries and medical conditions.

Principle

The urine alcohol test is a chemical analysis based on an alcohol-sensitive enzymatic reaction. If alcohol is present in the sample, it reacts with the chemicals on the reaction pad and causes a change in color. The Lab-Away urine alcohol test consists of a reaction pad. When in contact with alcohol solutions, the reaction pad rapidly changes color depending on the concentration of alcohol present. The pad uses solid-phase chemistry that utilizes a highly specific enzymatic reaction.

Reagents

Tetramethylbenzidine
Alcohol oxidase

Peroxidase
Other additives

Precautions

The Lab-Away Rapid Urine Alcohol Concentration Test is a visually interpreted test that uses color matching to obtain an estimate of the relative blood alcohol concentration.

Test materials that have been exposed to urine should be treated as potentially infectious.

Do not use the Lab-Away Rapid Urine Alcohol Concentration Test after the expiration date marked on the package.

Restrictions for compliance with appropriate preventive measures and local regulations.

Storage and Persistence

The Lab-Away Rapid Urine Alcohol Concentration Test should be stored at 2-8°C (36-46°F) in the sealed foil package. If the storage temperature exceeds 8°C (68°F), the test results may become less reliable. Before opening the package, the test should be placed at room temperature.

Materials provided

Individual test strip
Urine cup
Pipette for transferring test urine onto the strip

Materials not provided (but needed)

Stopwatch for timing prior to interpretation of results.

Instructions for use


Allow the strip to equilibrate to room temperature (15-27°C) before testing.

- 1) Refrain from placing anything in the mouth for fifteen minutes before beginning the test.
- 2) This includes soft drinks, tobacco products, coffee, breath mints and food, etc. Open the foil packet and remove the test strip.

3) Remove the test strip from the sealed package and place it on a clean and flat surface. The area of the strip on which the urine sample is applied should be facing up.

4) Using the pipette, transfer one drop (approximately 40 µL) from the urine cup onto the designated area of the test strip, taking care not to bubble.

5) Two minutes later, the results will be ready. You can interpret them by comparing the colors on the test strip with the corresponding value (in ppm) on the color card.

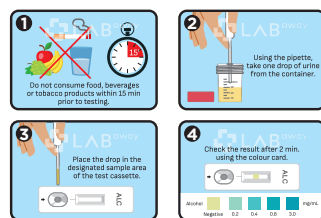

Color card



Interpretation of results

Positive: The response area changes from light blue to dark blue, falling on or between the corresponding color blocks on the color map.

Negative: No colour change appears on the reaction area. The color should correspond with the area on the color card that indicates a negative (-) result. This indicates that no alcohol was detected in the urine tested

Invalid: If the control zone on the test strip has a blue colour.

Limitations

1) The Lab-Away Rapid Urine Alcohol Concentration Test is very sensitive to the presence of alcohol. Occasionally alcohol vapors in the air are detected by the test. Alcohol fumes are present in many institutions and homes. Alcohol is an

ingredient in many household products such as disinfectants, deodorants, perfumes and glass cleaners. If you suspect the presence of alcohol fumes around you, test in a place known to be free of them.

2) Ingestion or use of medications from pharmacies or drugstores, as well as medications containing alcohol, may result in positive results even without knowing consumption of alcoholic products.

3) The concentration of alcohol in the human body slowly increases after alcohol ingestion. In most cases, the maximum concentration of alcohol in human urine occurs within 30 to 60 minutes after the last ingestion of alcohol. After the maximum occurs, the alcohol concentration in the body begins to decrease. How long it takes for the body's alcohol content to reach 0 depends on the amount consumed, gender, weight, alcohol tolerance, and a number of other factors.

Influencing substances

The following substances may affect the Lab-Away Rapid Urine Alcohol Test when using samples other than urine. The substances listed are usually not present in urine in sufficient quantity to affect the test.

A) Ingredients that enhance colour development

Peroxidases
Strong Oxidates

B) Constituents that inhibit colour development

Reducing agents: ascorbic acid, tannic acid, pyrogallol, mercaptans and tosylates, oxalic acid, picric acid.

Bilirubin

L-DOPA

L-Methyldopa

Methamphetamines

Bibliography

1. Diagnosis, Clinical Aspects and Biopsychosocial Causes, Substance Abuse Library, University of Pennsylvania, 1997.
2. Jones, A. W.: Inter- and intraindividual variations in the saliva/blood alcohol ratio during ethanol metabolism in man, Clin. Chem. 25, 1394-1398, 1979.
3. MaCall, L.E.L., Whiting, B., Moore, M.R. and Goldberg, A. : Correlation of ethanol concentrations in blood and saliva, Clin. Sci. 56, 283-286, 1979.

MANUFACTURER

Qingdao Hightop Biotech Co., Ltd.
369 Hedong Road, Hi-tech Industrial Development Zone,
Qingdao Shandong, 266112, China

DISTRIBUTER

„Gryphone Health Solutions
EOOD
30, Ekzarh Yosif str.
1000 Sofia
Bulgaria
Contact: +359 882 399 001
office@lab-away.com

Please check the instructions before using. Not reusable.
Store between 4 - 30°C



www.lab-away.com